

## Class-leading Controlled Heating

### HEATING

- ▮ 0.2 Seconds for 200 °C
- ▮ 0.5 Seconds for 250 °C
- ▮ Quick heating prevents PDL damage

### COOLING

- ▮ 2 Seconds for Cooling
- ▮ Swift cooling makes sure GP doesn't follow along the heating tip

*FAST PACK*



*FAST PACK*

## Internally Heated Tip

*Smart, Safe, Sophisticated*

- | Only the working end of the tip is heated- avoids need for insulated sleeves
- | The non-working end is insulated from within and stays cool- making it comfortable for patients and also saves battery
- | Safety cap where the tip's non-working end is plugged into the device- safer for you while changing tips




90° - 250° C

## Wide Temperature Range

Customize temperature and time as per obturation material used or any other intended use of the device

*FAST PACK*

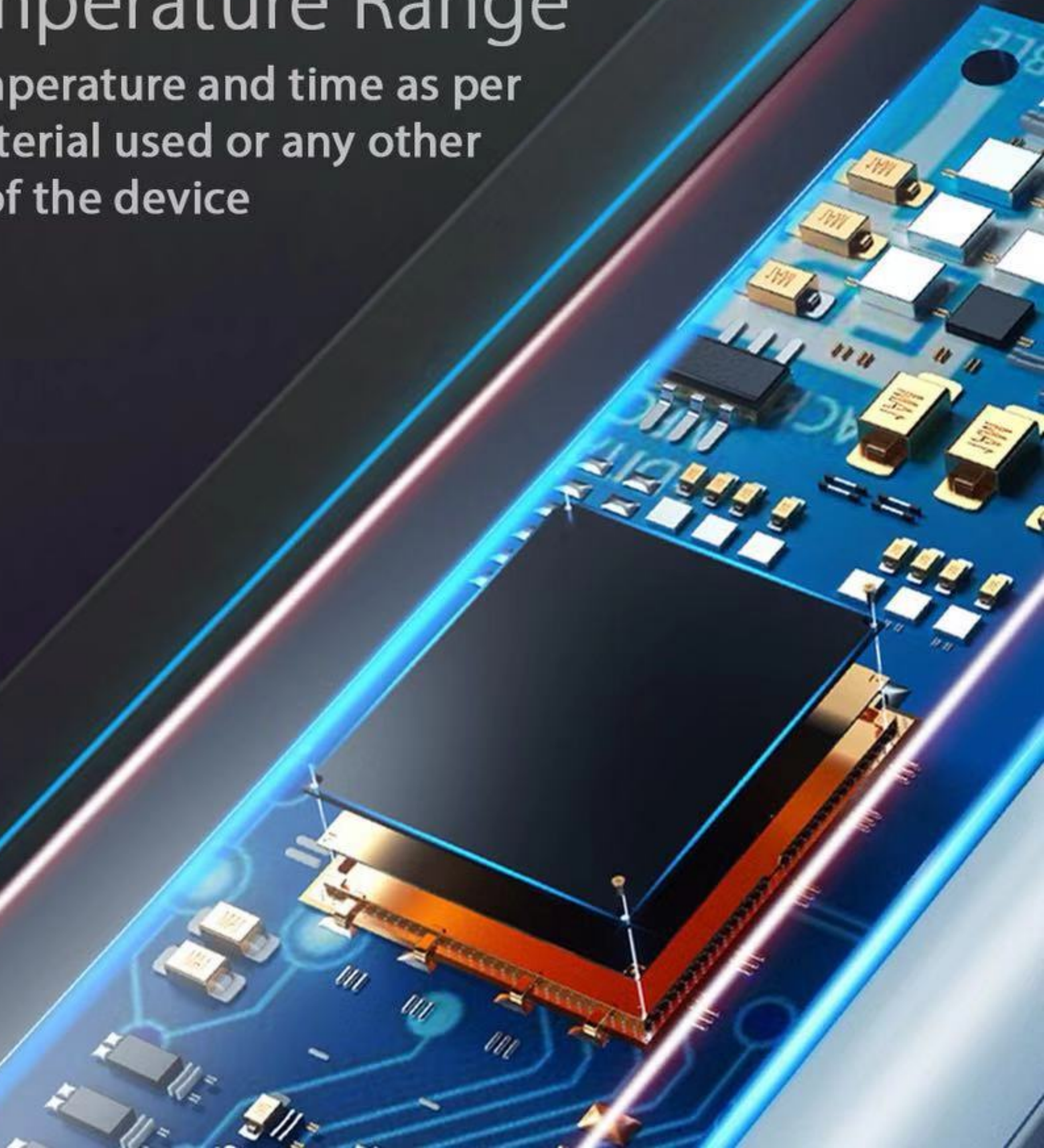

# Experience Ergonomics

- I 360° start/stop ring allows holding the device in an angle comfortable for you
- I Left-right handed operation
- I Comfortable grip
- I Light in weight

*FAST PACK*



*FAST PACK*

# Discover the Powerful Battery

- | Highly efficient Li-ion batteries which retain charge longer
- | 2 hours to complete charge
  - | 15 minute turbo charge to 30%
  - | Auto-switch to regular charging after that
- | 4 hours of continuous operation at 250°C
- | 50 cases with single charge
- | Option to direct plug-in to electric supply when zero battery

